DATA SHEET

Mass Flow Controllers & Meters

5800E Series

Elastomer Sealed, Analog, General Purpose Thermal Mass Flow Meters & Controllers for Gases



Model 5850E

The 5800E Series mass flow meters and mass flow controllers have gained broad acceptance as the standard for stability and reliability. These products have a wide flow measurement range and are suitable for a broad range of temperature and pressure conditions making them well suited for applications in chemical and petrochemical research, laboratory, analytical, fuel cell and life science among others.

Highlights of the 5800E Series mass flow products include: industry leading long-term stability, accuracy backed by superior metrology systems and methods using primary calibration systems directly traceable to international standards, and a range of analog I/O options.

The 5800E Series provides a highly configurable platform based on a simple modular architecture. The 5800E Series feature set was carefully selected to enable drop-in replacement of many brands of mass flow controllers. With the wide range of options and features available, the 5800E Series provides users with a single platform to support a broad range of applications.

Features	Benefits
Industry leading long-term sensor stability	Increased system uptime and reduced cost of ownership by reducing maintenance and eliminating periodic recipe adjustments and/or recalibrations
Superior valve technology	Minimum leak-by, wide turndown, and superior corrosion resistant materials reduces overall gas panel cost
Adaptable mechanical configurations	Easily retrofit to existing systems
Primary standard calibration systems	Ensures measurement accuracy is traceable to international standards
Simple modular design and reducing total cost of ownership	Easy-to-service elastomer sealed design provides for factory or field service maximizing uptime

View 5800 Product Page



Advanced Thermal Flow Measurement Sensor

Brooks' sensor technology combines:

- Excellent signal to noise performance for improved accuracy at low setpoints
- Superior long-term stability through enhanced sensor manufacturing and burn in process
- Isothermal packaging to reduce sensitivity to external temperature changes


Wide Flow Range

The 5800E Series covers an extremely broad range of flow rates. Model 5850 can have a full scale flow as low as 3 ccm. With a high turndown ratio of 50:1 repeatable gas flow can be measured or controlled down to 0.06 ccm! Model 5853 can monitor or control gas flows up to 1000 lpm.

Fast Response Performance

The electronics and superior mechanical configuration in the 5800E Series provide for fast, stable response characteristics.



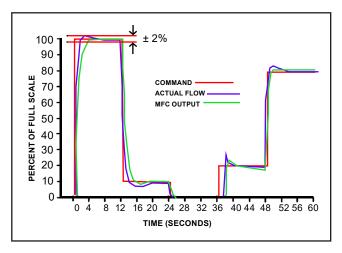

Flow Ranges and Pressure Ratings:

Mass Flow Controller	Mass Flow Meter		Ranges 4. Ratings	Pressure Unit psi/bar	PED Module H Category
Model	Model	Min. F.S.	Max. F.S.	Standard	
5850	5860	0.003	30 lpm	1500 psi/100 bar	SEP
5851	5861	10	100 lpm1	1500 psi/103 bar	SEP
5853	5863	100	1000 lpm	1000 psi/70 bar	1 for all 150 lb flanges 2 for all other connections

 $^{^{\}mathtt{1}}$ 200 lpm of H2 possible, 600 lpm of H2 possible with decreased accuracy

	5850/60	5851/61	5853/63	
PERFORMANCE				
Flow Accuracy ²	1% F.S.	1% F.S.	1% F.S.	
Control Range		50:1		
Repeatability & Reproducibility		0.25% of rate		
Linearity	Included in accuracy <3 seconds			
Response Time (Settling Time within ±2% F.S. for 0-100% command step)				
Zero Stability		$< \pm 0.2\%$ F.S. per year		
Temperature Coefficient	Zero: <0.075% of F.S. per °C. Span: <1.0% of F.S. shift from original calibration over 10-50°C (50-122°F) range			
Pressure Coefficient	±0.03% per psi (0-200 psi N2)			
Attitude Sensitivity	<0.5% F.S. maximum deviation from specified accuracy after re-zeroing			

² Accuracy including linearity at calibration conditions.

RATINGS

Operating Temperature Range 5-65°C (41-149°F)

Minimum Pressure Differential (Controllers) ³	5 psi/0.35 bar	10 psi/0.69 bar	Min.: 7.5 psi/0.52 <500 lpm Min.: 11.8 psi/0.81 >500 lpm
Maximum Pressure Differential (Controllers) ³	50 psi/3.45 bar		290 psi/20 bar
Leak Integrity (external)	1x10 ⁻⁹ atm. cc/sec He		

³ Differential pressures beyond the specified limits may be possible depending on process conditions.

N 4	1 1 1	N III	CAL


Valve Type Normally Closed, Normally Open, Meter

Primary Wetted Materials 316L Stainless Steel, High Alloy Stainless Steel, Viton® fluoroelastomers, Buna-N, Kalrez® and Teflon®


•				
ELECTRICAL	5850E	5851E	5853E	5860E/5861E/5863E
Electrical Connection	15-Pin Male Sub D-Type (DA-15P) or Card Edge connector (Controller models only)			
Input (Setpoint) Signal	0 to 5 Vdc (200 K ohms input resistance) N/A			
Output Signal	0 to 5 Vdc into 2000 ohms (or greater) load			
Reference Output Signal	5 Vdc ±0.01 Vdc, max. load 2 K ohm			
Power Requirements	N. C. Valve or N.O. Valve	10.5 watts	3.5 watts	1 05 watts


Reference Output Signal N. C. Valve or N.O. Valve with flow <2.5 slpm: 3.5 watts, +15 Vdc @ 35 mA, -15 Vdc @ 180 mA N.O. Valve with flow rate > 2.5 slpm: 10.5 watts, +15 Vdc @ 350 mA, -15 Vdc @ 350 mA	1 3				
with flow <2.5 slpm: +15 Vdc @ 35 mA, +15 Vdc @ 35 mA, 3.5 watts, +15 Vdc @ 35 mA, -15 Vdc @ 350 mA -15 Vdc	Reference Output Signal	5 Vdc ±0.01 Vdc, max. load 2 K ohm			
	Power Requirements	with flow < 2.5 slpm: 3.5 watts, +15 Vdc @ 35 mA, -15 Vdc @ 180 mA N.O. Valve with flow rate > 2.5 slpm: 10.5 watts, +15 Vdc @ 350 mA,	+15 Vdc@ 350 mA,	+15 Vdc @ 35 mA,	+15 Vdc @ 35 mA,

5850E, Card Edge Connector

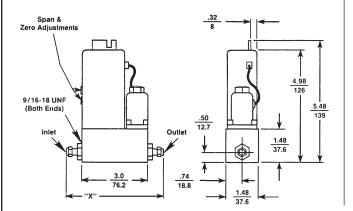

5850E, 15-Pin D Connector

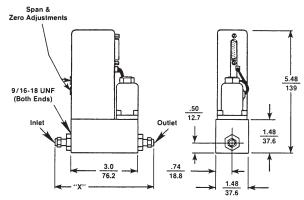


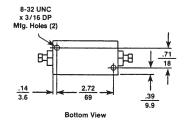

Connection Type	"X"Dim.
1/8"Compression Fitting	4.84/122.9
1/4" Compression Fitting	5.02/127.5
3/8" Compression Fitting	5.14/130.5
1/4" Tube VCO®	4.56/115.8
1/4" Tube VCR®	4.88/124.0

Inches/Millimeters

5860E, 15-Pin D Connector

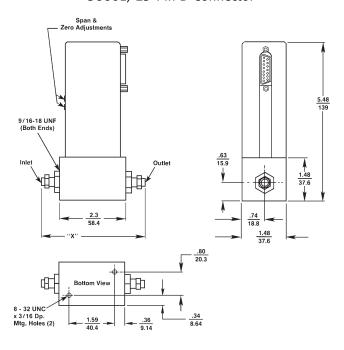


Connection Type	"X" Dim		
1/8" Compression Fitting	4.14/105.2		
1/4" Compression Fitting	4.32/109.7		
3/8" Compression Fitting	4.44/112.8		
1/4" Tube VCO®	3.86/98.0		
1/4" VCR®	4.18/106.2		
I I (A dillion - de o-			


Inches/Millimeters


Product Dimensions


5851E, Card Edge Connector

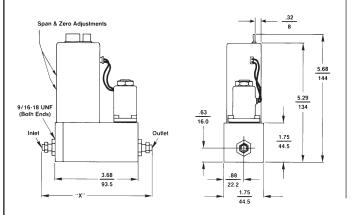

5851E, 15-Pin D Connector

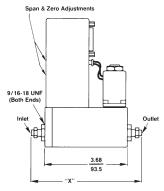

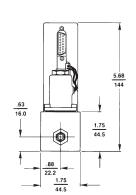


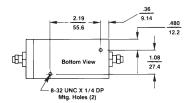



Connection Type	"X"Dim.
1/4"Compression Fitting	5.70/145
3/8" Compression Fitting	5.82/148
1/4" Tube VCO®	5.22/133
1/4" Tube VCR®	5.56/141
1/2" Tube VCO®	5.68/139
1/2" Tube VCR®	5.86/149

Inches/Millimeters

5861E, 15-Pin D Connector

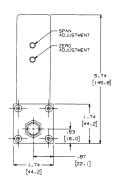

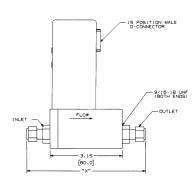


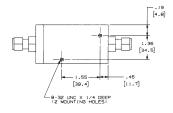


Connection Type	"X" Dim
1/4" Compression Fitting	5.17/131
1/4" Tube VCO®	4.71/120
1/4" Tube VCR®	5.03/128
3/8" Compression Fitting	5.29/134
3/8" or 1/2" Tube VCO®	5.15/131
3/8" or 1/2" Tube VCR®	5.33/135

Inches/Millimeters

5853E, Card Edge Connector With Fittings

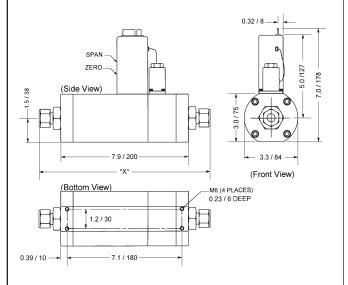


5853E, 15-Pin D Connector With Fittings

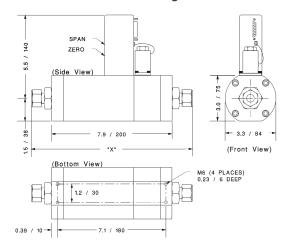


Connection Type	"X" Dimension	
	inch	mm
Female Pipe Thread	7.9	200
1/2" Tube Compression*	10.6	268
3/4" Tube Compression*	10.6	268
1" Tube Compression*	10.9	277
1/2" VCO** (100 slpm max.)	9.9	251
3/4" VCO*	10.2	258
1/2" VCR** (100 slpm max.)	10.1	257
1/2 VCR*	10.2	259

^{* 1-1/16 - 12} UNF Connection To Body

5853E, Card Edge With Flanges

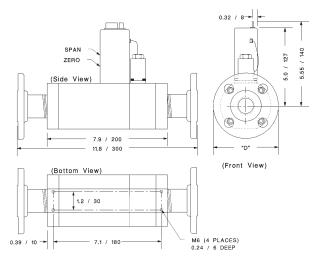

5853E, 15-Pin D Connector With Flanges

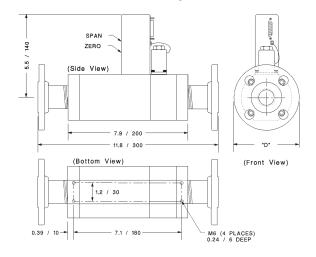


Flange	Flange OD, "D"		
Туре	inch	mm	
ANSI 1/2" 150 lb	3-1/2	89	
ANSI 1/2" 300 lb	3-3/4	95	
ANSI 1.0" 150 lb	4-1/4	108	
ANSI 1.0" 300lb	4-7/8	124	
ANSI 1.5" 150 lb	5	127	
ANSI 1.5" 300 lb	6-1/8	156	
ANSI 2.0" 150 lb	6	152	
ANSI 2.0" 300 lb	6-1/2	165	
DIN DN15 PN40	3.74	95	
DIN DN25 PN40	4.53	115	
DIN DN40 PN40	5.90	150	
DIN DN50 PN40	6.50	165	

Inches/Millimeters

^{** 9/16 - 18} UNF Connection To Body

5863E, 15-Pin D Connector With Fittings

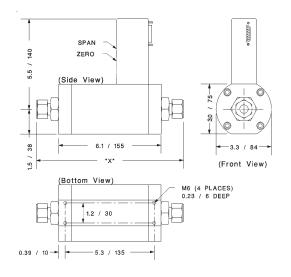


Connection Type	"X" Dimension	
	inch	mm
Female Pipe Thread	6.1	155
1/2" Tube Compression*	8.8	223
3/4" Tube Compression*	8.8	223
1" Tube Compression*	9.1	232
1/2" VCO ** (100 slpm max.)	8.1	206
3/4" VCO *	8.4	213
1/2" VCR ** (100 slpm max.)	8.3	211
1/2" VCR*	8.4	214

* 1-1/16 - 12 UNF Connection To Body

** 9/16 - 18 9/16 - 18 UNF Connection To Body

5863E, 15-Pin D Connector With Flanges

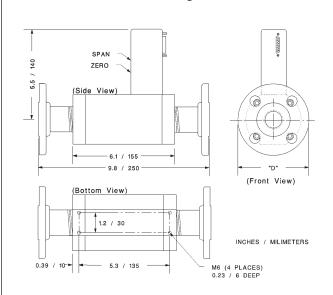

Flange	Flange OD, "D"		
Туре	inch	mm	
ANSI 1/2" 150 lb	3-1/2	89	
ANSI 1/2" 300 lb	3-3/4	95	
ANSI 1.0" 150 lb	4-1/4	108	
ANSI 1.0" 300lb	4-7/8	124	
ANSI 1.5" 150 lb	5	127	
ANSI1.5" 300 lb	6-1/8	156	
ANSI 2.0" 150 lb	6	152	
ANSI 2.0" 300 lb	6-1/2	165	
DIN DN15 PN40	3.74	95	
DIN DN25 PN40	4.53	115	
DIN DN40 PN40	5.90	150	
DIN DN50 PN40	6.50	165	

Inches/Millimeters

Certifications

Certifications

Mark	Agency	Certification	Applicable Standard	Details
ϵ	CE	EMC Directive 2004/108/EC	EN:61326-1:2006	PASS

Service and Support

Brooks is committed to assuring all of our customers receive the ideal flow solution for their application, along with outstanding service and support to back it up. We operate first class repair facilities located around the world to provide rapid response and support. Each location utilizes primary standard calibration equipment to ensure accuracy and reliability for repairs and recalibration and is certified by our local Weights and Measures Authorities and traceable to the relevant International Standards.

Visit www.BrooksInstrument.com to locate the service location nearest to you.

START-UP SERVICE AND IN-SITU CALIBRATION

Brooks Instrument can provide start-up service prior to operation when required. For some process applications, where ISO-9001 Quality Certification is important, it is mandatory to verify and/or (re)calibrate the products periodically. In many cases this service can be provided under in-situ conditions, and the results will be traceable to the relevant international quality standards.

SEMINARS AND TRAINING

DS-TMF-5800E-Series-MFC-eng/2023-1

Brooks Instrument can provide customer seminars and dedicated training to engineers, end users, and maintenance persons. Please contact your nearest sales representative for more details. Due to Brooks Instrument's commitment to continuous improvement of our products, all specifications are subject to change without notice.

TRADEMARKS Brooks Brooks Instrument, LLC All other trademarks are the property of their respective owners.





Global Headquarters

Brooks Instrument 407 West Vine Street Hatfield, PA 19440-0903 USA Toll-Free (USA): 888-554-FLOW T: 215-362-3500

BrooksAM@BrooksInstrument.com

A list of all Brooks Instrument locations and contact details can be found at www.BrooksInstrument.com

©Copyright 2023 Brooks Instrument, LLC All rights reserved. Printed in U.S.A.

